

# EDMUS 5.7.1 – INSTALLATION INSTRUCTIONS

For a **network installation**, see the accompanying document:  
"EDMUS 5.7.1 - Network Installation Instructions.pdf".

To be able to make a NEW installation of EDMUS,  
**you must have administrator rights on the computer** (for the duration of the installation).

(Note: the illustrations were made for EDMUS 5.7, they will mention 5.7.1 in the new version)

## BEFORE THE INSTALLATION

### 1. Important: first make a backup copy of your EDMUS database


If you already have an EDMUS installation,  
**make sure you have a recent backup copy of your EDMUS database**  
(the file containing your patients' data) before proceeding with the installation.

The database file is usually named "Edmus.df1", but may have been renamed.

#### If you have an EDMUS 5 or higher database


- *Make a backup copy directly from EDMUS*

Choose "Back up Database" from the  
"Manager" menu to **make a backup copy  
of your database file** at the location of  
your choice.

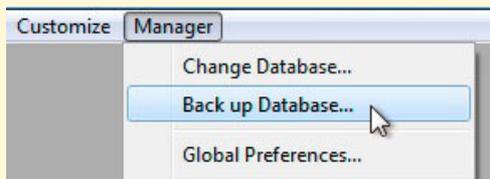

- *Make a backup copy manually*

To find out where your EDMUS datafile is  
located, choose "About EDMUS" from the  
"File" menu.:

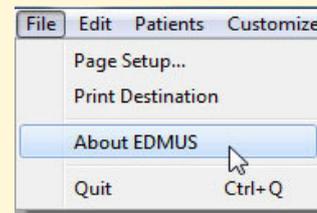


In the window that opens, the path to the  
database file is shown:

```
Database C:\USERS\PUBLIC\EDMUS SOFTWARE  
V5\EDMUS\data\Edmus.df1
```

**Make a backup copy of this file** before  
continuing with the installation.


Important: if there is also a file with the  
same name as the database file and a .DF2  
extension, copy it as well: it is part of the  
database.

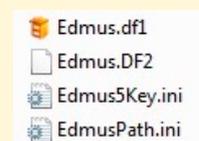


### If you have an EDMUS 4 or 4.1 database

- o *Make a backup copy manually*

To find out where your EDMUS datafile is located, choose "Change Database" from the "Manager" menu:

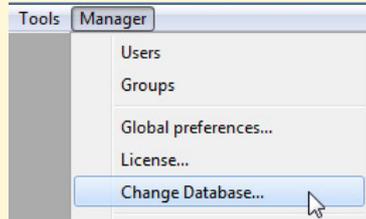

In the window that opens, the path to the database file is shown:

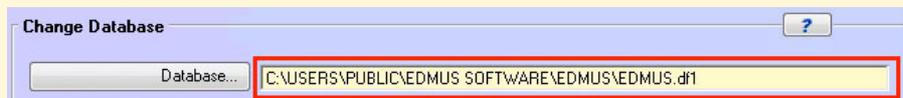

- o **Make a backup copy of this file** before continuing with the installation.

### Do NOT continue until you have a backup copy of your EDMUS database file.

Keep this copy in a safe place, at **another** location (external hard disk, USB stick etc.) than the original database.

If you cannot find the database file to back up, please contact technical support ([support@edmus.org](mailto:support@edmus.org)).

## 2. If you haven't already, download EDMUS 5.7.1 from the EDMUS web site

[http://www.edmus.org/en/soft/edmus\\_dl.php](http://www.edmus.org/en/soft/edmus_dl.php)




The downloaded "EDMUS 5.7.1 - Installation.zip" file is a zipped folder, you should unzip it:

- o *Windows*: right-click on the file, then choose "Extract" or "Extract All" (the wording may vary with your version of Windows).
- o *Mac*: double-click the file.

The "EDMUS 5.7.1 - Installation" folder will be created next to the .zip file.

The "EDMUS 5.7.1 - Installation" folder contains:

- o the installation program itself ("EDMUS 5.7.1 installer"),
- o and installation instructions in English and French.

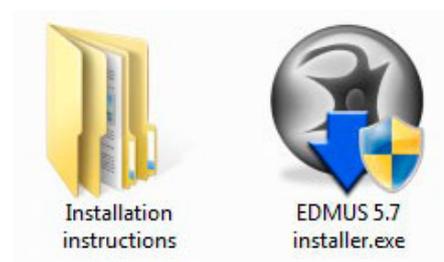

## THE INSTALLATION


### 3. Launch the installation program

Double-click the installation program ("EDMUS 5.7.1 installer") to start the installation of EDMUS.



#### A note on security settings

When launching the EDMUS installer program for the first time, depending on the security settings of your computer, you may encounter a message telling you that the installer program comes from "an unknown publisher" or "an unidentified developer":

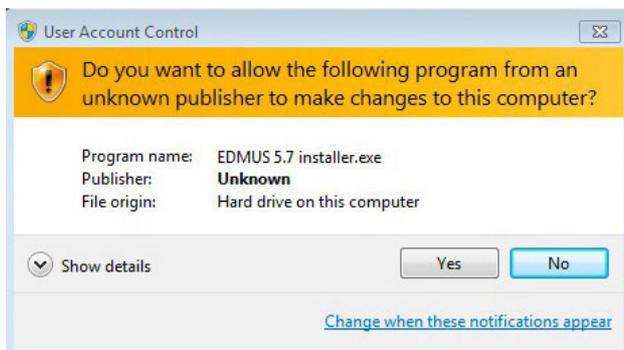

In this case:

*On Windows*

Click "Yes" to "allow the following program from an unknown publisher to make changes to this computer".

*On Mac OS X*

Right-click on the installer and choose "Open", then click "Open" in the dialog window that appears.

## 4. Proceed with the installation

---

When installing, three situations may occur:

### CASE 1: You are installing EDMUS for the first time on this computer

No option will be proposed. Accept the terms of the license, then let the installation proceed.


#### What's going to happen

EDMUS 5.7.1 will be installed, with an empty database (no patient).

### CASE 2: You already have a previous version of EDMUS 5 (5.0 to 5.7) on this computer

Accept the terms of the license, then select the "UPGRADE" option (it will be selected by default).

Then let the installation proceed.

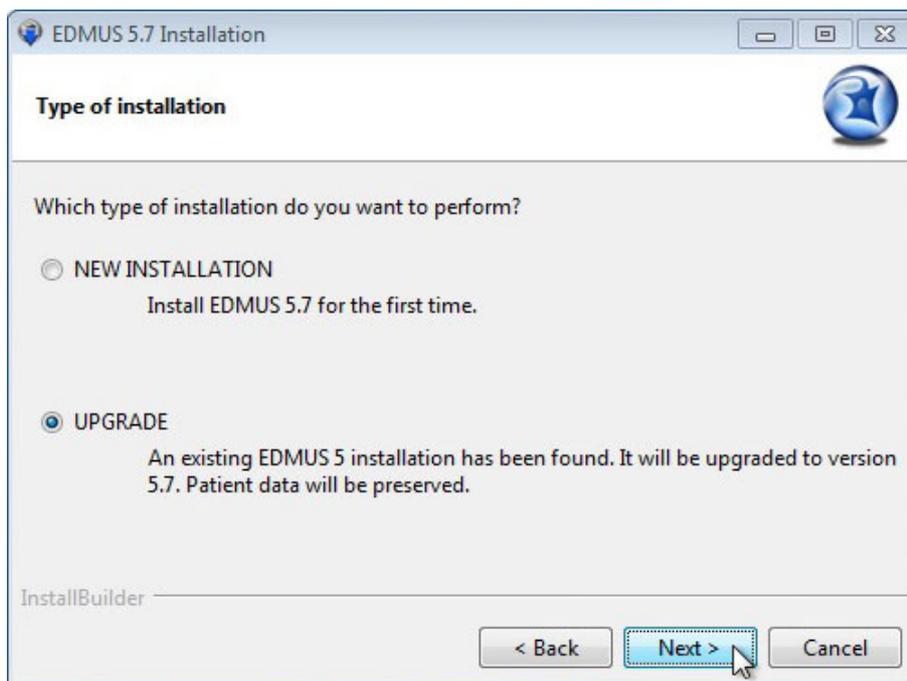


#### What's going to happen

EDMUS 5.7.1 will be installed, replacing EDMUS 5 (5.0 to 5.7). Your EDMUS datafile (with your patients' data) will be untouched. It will be updated at first launch of EDMUS 5.7.1 (see step 6 below).

### CASE 3: You already have EDMUS 4.0 or 4.1 on this computer

Accept the terms of the license, then select the "UPGRADE" option (it will be selected by default).

Then let the installation proceed.

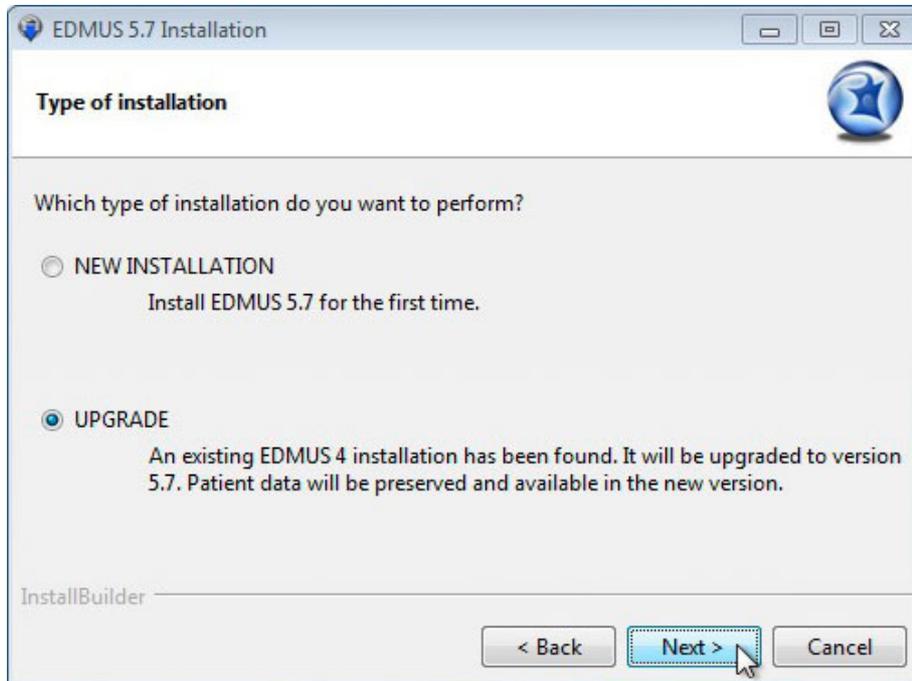

#### ***What's going to happen***

EDMUS 5.7.1 will be installed, without replacing EDMUS 4.0 or 4.1. Your EDMUS 4 database (with your patients' data) will be copied into the EDMUS 5.7.1 installation.

This copy of the datafile will be updated upon the first launch of EDMUS 5.7.1 (see step 6 below).

***At the end of the installation***, the installer program will recommend that you uninstall EDMUS 4:

- In most cases, leave the box checked to uninstall EDMUS 4, and click "Finish".
- However, **if you have no recent backup copy of your EDMUS 4 database**, uncheck the "Uninstall EDMUS 4" box. In this case, make the backup copy of the database you should have made at step 1 above.



You will then have to uninstall EDMUS 4 yourself:

- On Windows: Click the Start button, then "Control Panel", then "Add or Remove Programs" (the wording may vary slightly with different versions of Windows). Select "EDMUS 4", then click "Change/Remove".
- On Mac OS X: Drop the "/Applications/EDMUS Software" folder into the Trash (attention, **NOT** the "EDMUS Software v5" folder!). Empty the Trash.

## AFTER THE INSTALLATION

### 5. Once the installation is done, launch EDMUS 5.7.1 for the first time

#### o On Windows

Launch EDMUS 5.7.1 by double-clicking the "EDMUS 5" shortcut that was placed on the Desktop.



#### o On Mac

Launch EDMUS 5.7.1 by clicking the "EDMUS 5" shortcut that was placed in the Dock.


### 6. Update your older EDMUS database if any

If there was an older version of EDMUS on the computer before the installation of EDMUS 5.7.1, the database file (containing your patients' data) needs to be updated.

Two cases may occur depending on the version of EDMUS you had, corresponding to CASES 2 and 3 above:

#### CASE 2: There was an older version of EDMUS 5 (5.0 to 5.7)

At first launch of EDMUS 5.7.1, you will be asked if you wish to **update the database to the new EDMUS 5.7 format**, click Yes:

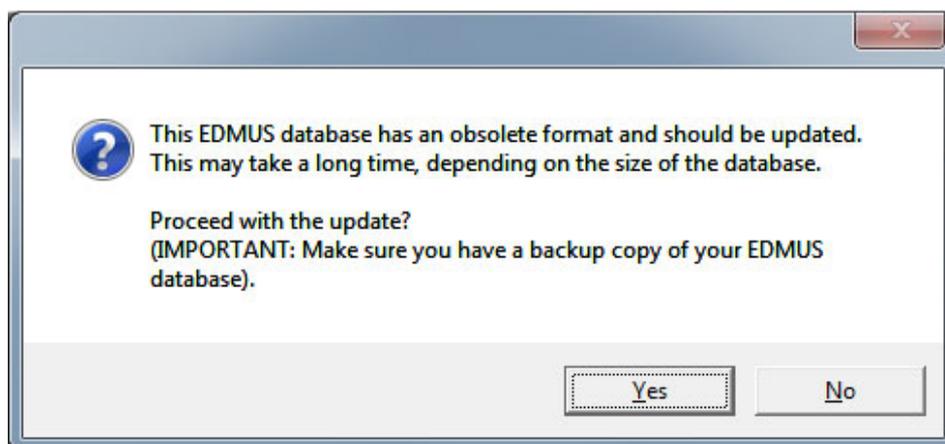

#### **What's going to happen**

EDMUS 5.7.1 will automatically update the database structure.


Attention, this operation can take a long time, up to one hour or more for a very large database (several thousand patients). This is very dependent on how powerful your computer is.

#### *IMPORTANT: Network installation*

If you are upgrading a network EDMUS installation, you should upgrade a LOCAL copy of the database, before placing it back onto the server. For the full procedure, see the accompanying document: "EDMUS 5.7.1 - Network Installation Instructions.pdf".

### CASE 3: There was an EDMUS 4 database (4.0 or 4.1)

- At first launch of EDMUS 5.7.1, the EDMUS 4 database will automatically be **converted to Unicode** (if this step does not proceed automatically, choose the "Full conversion" option when asked):

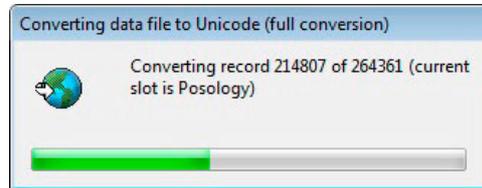

- Then (exactly as in CASE 1 above), you will be asked if you wish to **update the database to the new EDMUS 5.7 format**, click Yes.

#### **IMPORTANT: Network installation**

If you are upgrading a network EDMUS installation, you should upgrade a LOCAL copy of the database, before placing it back onto the server. For the full procedure, see the accompanying document: "EDMUS 5.7.1 - Network Installation Instructions.pdf".

## 7. User name and Password

Au lancement (après la mise à jour de la base le cas échéant), EDMUS demandera *votre nom d'utilisateur et votre mot de passe*. Par défaut, il s'agit de ADMIN et ADMIN (ne pas tenir compte des majuscules ou minuscules), sauf si vous les aviez modifiés dans la base EDMUS antérieure que vous venez de mettre à jour.



## 8. At first launch: importing new Thesaurus and Help files

### ONLY WHEN UPGRADING AN OLDER EDMUS 4 OR 5 VERSION (CASE 2 or 3 above)

At first launch, EDMUS will **automatically import new Thesaurus and Help files**: this will only occur only this once, and will last a few minutes at most. The personal items you may have added to the Thesauruses will be kept.



Note: if, when launching EDMUS later, this importation step occurs again, you should **DELETE** the following files manually which may have remained in the folder where EDMUS was installed (EDMUS Software v5/EDMUSdata/):

```
Countries.txt, Departements.txt
EDMUSV5_EX.ETH, EDMUSV5_ST.STH, EDMUSV5_TH.THE, EDMUSV5_TT.TTH
EDMUSV5.HLP, EDMUSV5.SCL
```

## 9. Activate your EDMUS 5 license

Once EDMUS is installed and the database (if any) is updated,  
**you should activate your EDMUS 5 license.**

Without a license, you will not be able to create more than 3 patients in the database,  
 and EDMUS will stop functioning after a given date.

If you don't have an EDMUS 5 license, or you have lost it, please contact the technical support ([support@edmus.org](mailto:support@edmus.org)). The license is provided as a .zip file.

To activate your EDMUS 5 license:

- o Unzip the file, you will get one or more file(s) with an ".ini" extension (for instance "B1286.ini"), each corresponding to an EDMUS license.
- o Launch EDMUS 5, and select "About EDMUS" from the "File" menu.
- o In the window that opens, click the "Activate license" button.
- o Go select the license file.
- o You will then see that the "EDMUS serial number", "Licensed to" and "Database ID" items have been updated with your data.




## 10. EDMUS Language & Configuration

### o **Language**

At first launch, for a new EDMUS installation, the program will by default be in English.

If you wish EDMUS to use another language:

- Choose "User Preferences" from the "Customize" menu.

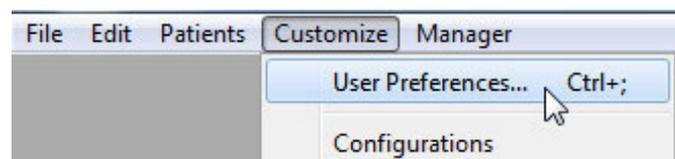

- ▶ Then, in the "User Preferences" window that appeared, click "Modify", select the language to use, and click OK.




## ○ **Configuration**

A configuration is a **presentation format**, a way of presenting data in a more or less detailed form or for a specific use.

At first launch, for a new EDMUS installation, the program will be by default in the "Light" configuration, which only displays the main fields of the database.

You can change the configuration:

- ▶ Choose "User Preferences" from the "Customize" menu.
- ▶ Then, in the "User Preferences" window that appeared, click "Modify", select the configuration to use, and click "OK".

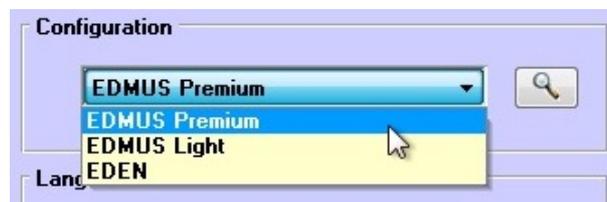

The following configurations are available:

- ▶ **Premium**: the most complete (displays all fields of the database),
- ▶ **Light**: simplified (only displays the main fields of the database),
- ▶ **EDEN**: specifically tailored for NMOSED patients.

In every case, the same database is used, the difference is only one of presentation.

### **EDMUS 5.7.1 - LAST VERSION OF EDMUS**

**Version 5.7.1, released on 2019-05-06, will be the last version of EDMUS.**

Development of the software has stopped and no further version will be provided.

Technical support remains available.

Important: Please note that, although EDMUS 5.7.1 is compatible with Windows 10, it is not certified on this platform and may not function optimally.